



LDA<sup>®</sup>, Automated Dynamic Laximetry

### Knee ligament analysis

- in tibial translation
- in tibial rotation







The company GENOUROB specializes in design, production and marketing of medical devices for the evaluation of the state and performance of the ligamentous structures of the knee.

We invite you with this document to discover the method LDA®, Automated Dynamic Laximetry, with its devices and their patented innovations.

We thank you for your interest and we remain at your disposal for any further information.

**Stéphane Nouveau**  
President and CEO

# Automated Dynamic Laximetry measurements / LDA<sup>®</sup> in automated tibial TRANSLATION



## ■ The LDA<sup>®</sup> method a fully automated Lachman's test

The innovation "Automated Dynamic Laximetry" (*Laxity Measurement*) offers an objective and reproducible method to perform Lachman's test, allowing an evaluation with accurate figures of ligament laxity of the knee at an automated translation.

## ■ A simple and fast test providing precise results

Once patient's data registered, the fixing parameters of ankle and femur stored, the sensor is placed on the ATT (Anterior Tibial Tuberosity), which measures the forward translation of the tibia, provoked by an automated push on the calf (proximal part).

By comparison of the measurements of both knees, the LDA<sup>®</sup> software displays the laxity curves, calculates the differential of the displacements and of the slopes.

## ■ NEW: inclusion of the slopes parameter in the functional analysis

So far only measuring the laxity differential, without consideration of the gradients of the curves (inverse of the stiffness) did not allow a comprehensive analysis of the ligament status (Bercovy and Weber\*). The LDA<sup>®</sup>, during translation pushes from 1 to 200 N (300 N maximum), registers more than 50 values of displacement, thus establishing accurate elongation curves with calculated slopes, giving objective assessment of the state of resistance of the Anterior Cruciate Ligament (H. Robert\*\*).

The GNRB detects incomplete and complete ACL ruptures and allows a functional analysis of the ligament.

## ■ NEW: induced medial rotation

The detection of incomplete ACL injuries is optimized by registering the medial rotation coupled to the anterior tibial translation. In this, a variation of rotation is measured, validated, thus indicating the risk of a partial ACL injury (P. Christel\*\*\*). This innovation characterizes the GNRB Rotab.

## ■ PCL injuries

Lesions of the Posterior Cruciate Ligament are likewise detectable with the LCP module (option), providing an automated posterior translation.


## ■ GNRB, more than just simple Laximetry!

Many technical innovations such as a biofeedback system avoiding false negatives, characterize the GNRB to optimize the accuracy and the reproducibility of the tests.

\* M. Bercovy & al. RCO; 1995; 81, 114-127,

\*\* H. Robert, & al. OTSR.2009; 95, 171-176

\*\*\* P. Christel & al. J Bone Joint Surg Br 2012; 94-B: 68-74

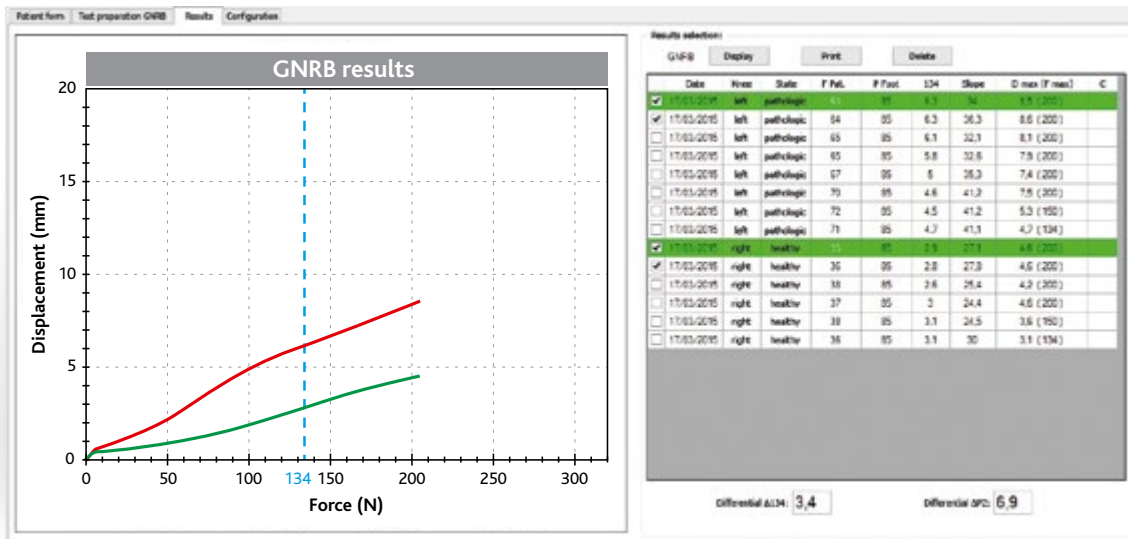

## Automated tibial translation parameters GNRB - GNRB Rotab (H. Robert, Isakos, Rio 2011)

Differential of displacements at 134 N ( $\Delta 134$ in mm)	Ligament state
$\Delta 134 > 3$	Complete lesion (according to $\Delta P2$ )
$1 < \Delta 134 < 3$	Partial lesion (according to $\Delta P2$ )
$\Delta 134 < 1$	No lesion


  

Differential of slopes ( $\Delta P2$ in $\mu\text{m}/\text{N}$ )	Risk of functional instability
$\Delta P2 > 10$	High
$5 < \Delta P2 < 10$	Medium
$\Delta P2 < 5$	Low

## At induced medial rotation / GNRB Rotab


Differential of rotations ( $\Delta^\circ$ in degree)	Risk of lesion of PL bundle
$\Delta^\circ > 3$	High
$\Delta^\circ < 3$	Low

**Knee**  
— Healthy    — Pathologic    — Operated

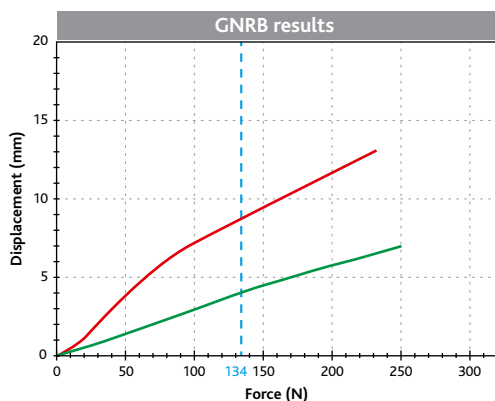


### Preoperative control

$\Delta 134 = 5,2$ ,  $\Delta P2 = 15$ ,  $\Delta^\circ = 2$   
 Objectification of a clinical Lachman's test (GNRB Rotab)  
 Complete rupture (no resistance to the push; high inclination...)

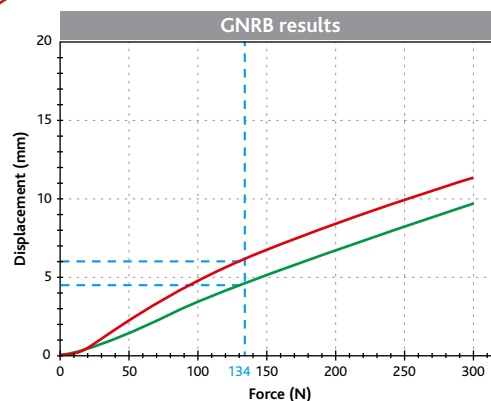

Similar  $\Delta 134 = 1,5$  mm but different slopes of curves  
 ↓  
**DIFFERENT FUNCTIONAL ANALYSIS**

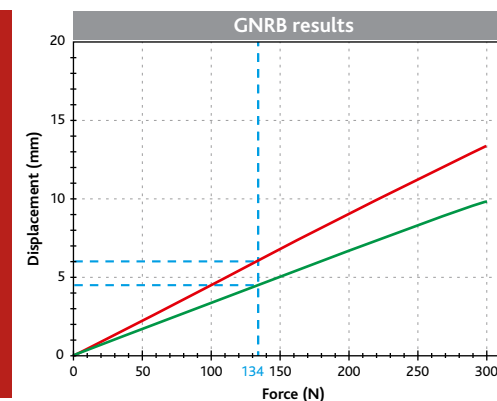

For identical  $\Delta 134$ , only the inclusion of the differential of the slopes permits an efficient functional analysis of the ligament state. This one is different according to the differential of the slopes  $\Delta P2$  (parallelism or divergence of the curves) associated with the differential of the displacements  $\Delta 134$ .

To the left: no risk of functional instability

To the right: risk of functional instability

# Automated Dynamic Laximetry measurements / LDA<sup>®</sup> in controlled tibial ROTATION

## ■ An additional evaluation of the peripheral knee ligament structures

The innovation of the LDA<sup>®</sup> measurement of controlled tibial rotation is to permit an objective and rigorous evaluation of the rotational laxities.

The LDA<sup>®</sup> method allows a new approach to lesional damages of the peripheral ligament structures that may be involved in the rotational instability.

## ■ A simple and fast method

Patient data is stored and the fixation parameters of the ankle and the femur are saved, a motor-torque of 1 to 8 Nm is applied to the tibia-ankle-foot block with a registration of the provoked rotation, as well medial as lateral.


The very strict and tight fixation of the tibia ensures a clear measurement of the tibial rotation without reading errors by parasite movements of the foot joints.

By comparing the measurements, the LDA<sup>®</sup> software displays the ligament elongation curves and the differential of rotation of both knees.



## ■ Automated measurements, reproducible and accurate!

Many technical innovations in the laximetry, such as the Biofeedback system detecting the hamstring muscle contractions (risk of false negatives), also characterize the ROTAM.

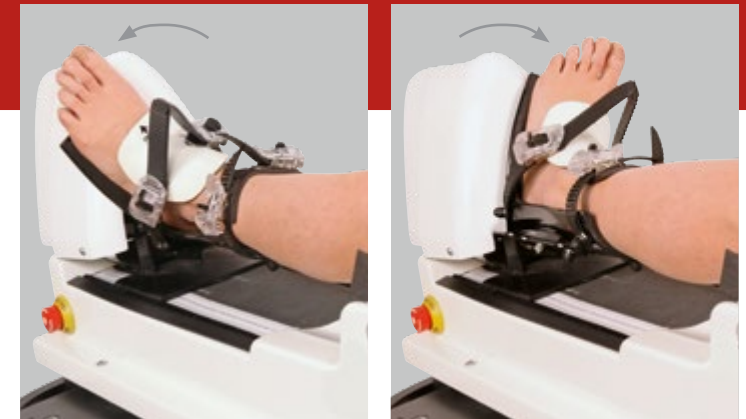

## ■ Objective preoperative assessments

Clinical or by MRI, lesions of peripheral ligament structures are sometimes difficult to verify and to quantify accurately.

The test in controlled tibial rotation stresses these structures to detect the peripheral lesions.

**In medial rotation:** an important differential of rotation indicates an affected ACL and anterior-lateral tibial structure (possible choice of an additional extra-articular surgery in association with ACL).

**In lateral rotation:** an important differential of rotation indicates an affected ACL and posterior-lateral tibial structure (possible choice of an additional extra-articular surgery in association with ACL).

## ■ Controlling postoperative results


Postoperative, the LDA<sup>®</sup> test of controlled tibial rotation demonstrates the efficacy of the "rotational brakes" and the quality of repairs at the ACL and at the peripheral structures.

## ■ The complementary partners ROTAM and GNRB

The ROTAM is the first motorized arthrometer dedicated to the objective dynamic evaluation of rotational laxity. It is the ideal complement to the GNRB, the automated laximeter to measure sagittal laxity.





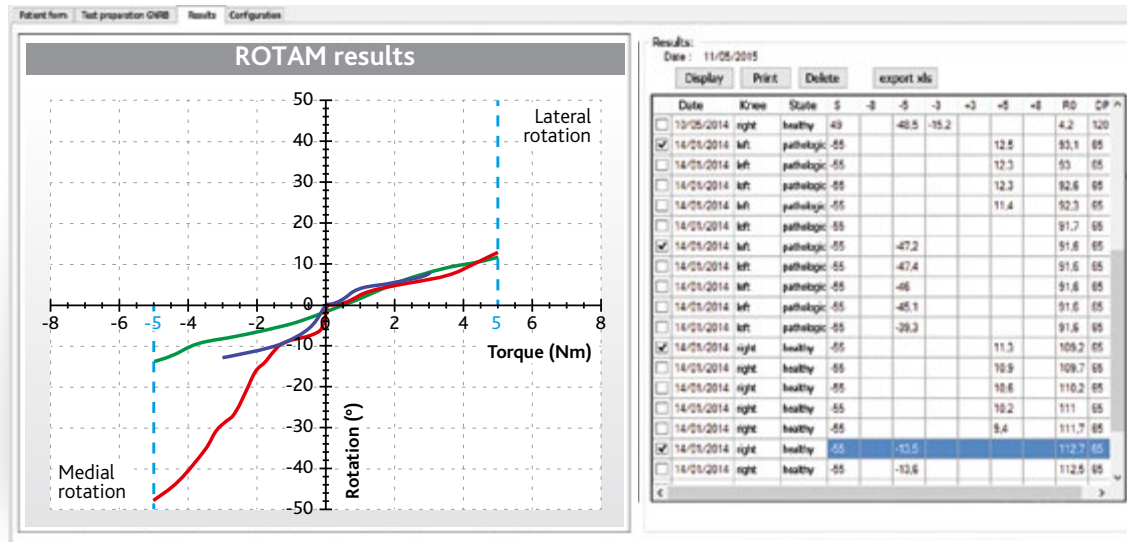


Test ROTAM at 5 Nm in controlled rotation	Differential of rotation $\Delta^{\circ}r > 6^{\circ}$	Differential of rotation $\Delta^{\circ}r < 6^{\circ}$
Analysis	Severe affection of the anterior-lateral peripheral structures	No severe affection of the anterior-lateral peripheral structures
Surgery (plastic) extra-articular	Recommended	Not recommended

(O. Lorbach, KSSTA, 2011 / H. Robert, ESSKA-SFA, 2014)




## Some examples:

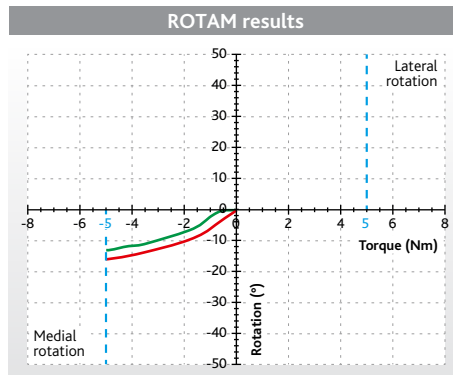


**Preoperative**

$\Delta^{\circ}r < 6$   
Normal state  
of peripheral ligament structures

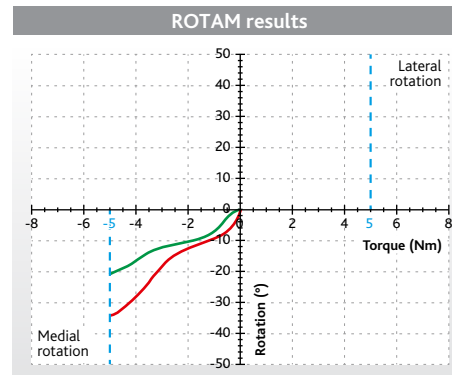


$\Delta^{\circ}r > 6$   
Divergent curves  
Pathological state of ligaments  
Extra-articular surgery recommended

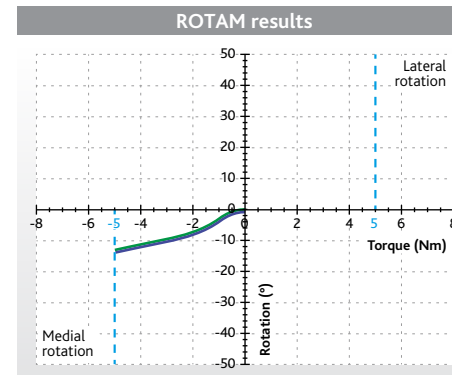

**Postoperative**

$\Delta^{\circ}r = 0$   
Successful reconstruction  
Parallel curves  
Good ligament resistance

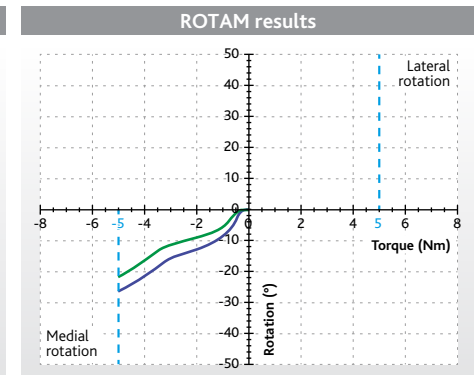

$\Delta^{\circ}r > 6$   
Parallel curves,  
but high risk of relapse

# The devices for the Automated Dynamic Laximetry



LCP  
GnrB



rotam



ROTAB  
GnrB

The range of our devices also includes the REHAB and the IP3, all of them can be combined and connected to the autonomous LDA® Station.



# Technical details and innovations



## Progressive tibial translation

Adjustable (1 to 300 N) or pre-selected (134 N, 150 N,... ) the progressive pushing force adapts to the individual need and guarantees a comfortable test for the patient.

## High precision laximetry

The sensor positioned on the anterior tibial tuberosity records 1/10 of a mm, and at each push, more than 50 tibial displacement values! An unmatched dynamic accuracy!



## Reproducible laximetry comparison

The patellar shell cap fixates the femur (patella against the femoral trochlea) with an individualized force, registered for each patient. This force applied identically to both knees warrants the comparison of the measurements and the reproducibility of the tests.

## Induced tibial rotation GNRB Rotab

The integrated electronic goniometer measures the tibial rotation induced by the anterior translation of the tibia. The distance of the foot positioning on GNRB and ROTAM devices is saved for a better reproducibility of tests.



## Controlled medial and lateral tibial rotation ROTAM

A torque force (in Nm) is progressively applied with the help of a motor to the unit tibia-ankle-foot. The device measures the provoked medial and lateral rotation to the degree and compares them with the values of the opposite knee.

## Accessories and options



### Station LDA®

With centralized control for several devices at choice and for the LDA® couch.



### Assistant LDA®

For PC, Printer and Device.



### Specific electrical couch for LDA®

With special guide plate for GNRB and ROTAM devices and a removable leg section for configuration in the examination couch.

## Comparison chart and characteristics

<b>Technical details</b>			
Type of Measures	Translation	Translation + medial rotation	Medial + lateral rotation
Accuracy	0,1 mm	0,1 mm / 1 Degree	1 Degree
Method of measurement in 50 points	✓	✓	✓
Manual selection of force / torque	✓ from 1 to 300 N	✓ from 1 to 300 N	✓ from 1 to 8 Nm
Automatic preselection of force / torque	✓ 134 - 150 - 200 - 250 N	✓ 134 - 150 - 200 - 250 N	✓ 3 - 5 - 8 Nm
Automatic Repetition of reading	✓ 1 to 3	✓ 1 to 3	✓
Detection and registration of muscular contractions (Biofeedback)	Option	Option	Option
Automatic registration of femur fixation force	✓	✓	✓
Registration of leg position	✓	✓	✓
Registration of patient data file	✓	✓	✓
Specific LDA® software	✓	✓	✓
Maintenance software integrated	✓	✓	✓
Dimensions and weights	845 x 270 mm (15 kg)	845 x 270 mm (17 kg)	845 x 270 mm (20 kg)
<b>Results analysis</b>			
Chart of measured values of displacements and rotations	✓	✓	✓
Display of Laximetry curves	✓	✓	✓
Calculation of rotation and translation differential between both knees	✓ Translation	✓ Translation and rotation	✓ Medial and lateral rotation
Calculation of inclination gradient and differential of the curves	✓	✓	✓
Printing and data transfer by network	✓	✓	✓
<b>Various options</b>			
Specific electrical couch for LDA®	✓	Required	Required
PC / Notebook - Printer	✓	✓	✓
Assistant LDA®	✓	✓	✓
<b>Quality Standards</b>			
Warranty	ISO 9001: 2008, ISO 13485 2 years	ISO 9001: 2008, ISO 13485 2 years	ISO 9001: 2008, ISO 13485 2 years

#### Quality Certificates

- NF IN ISO 13485 (2012)
- ISO 9001 (2008)
- ISO 13485 (2003)

#### Patents

- French patents (INPI): FR 0608725 and FR 0608726
- European patent: EP 078209.0-1526
- USA patent: Nr. 13/502790



Build. 60  
Rue du Chef de Bataillon Henri G eret  
53000 Laval - France

+33 (0)2 43 90 43 01  
contact@genourob.com

[www.genourob.com](http://www.genourob.com)